




**Digital Eye Chart**  
**3103 Sandstone Drive**  
**Idaho Falls, Idaho 83404**



**An affordable digital eye chart solution.**  
[www.digitaleyechart.com](http://www.digitaleyechart.com)  
 (208) 520 4292  
[sales@digitaleyechart.com](mailto:sales@digitaleyechart.com)



**An affordable digital eye chart solution.**  
[www.digitaleyechart.com](http://www.digitaleyechart.com)  
 (208) 520 4292  
[sales@DigitaleyeChart.com](mailto:sales@DigitaleyeChart.com)

The image on the front cover is the first of three installations of the Digital Eye Chart at The Vision Center of Jackson, Wyoming (two exam lanes and a pre-exam room). Enhance your patient's experience with a high tech alternative to projectors or paper charts.

**Chart Type:**  
Select Red/Green, Fixation, Four dot,

**Control:**  
Show all, Blank all, and darken background

**Choose Symbol Type:**  
Letters  
Numbers  
Symbols  
Tumbling E's  
Rolling C's

**Select Acuity:**  
Use menu to select an acuity or a set of 3 acuity rows to display

**Single Letter:**  
Click on any letter to show only one letter.

**Toggle Row:**  
Turn a row off and on.

**One Row Only:**  
Display this row only


**Randomize:**  
Randomize the chart with one click.

**One Click Blank:**  
Darken the chart to reduce illumination in the exam room. Return to chart by clicking again.



The doctor's laptop display has a control window with simple buttons and menu controls for controlling the chart the patient views on a wall mounted monitor. The doctor can easily check the patient's response without turning or interpreting the mirrored image on a projector screen. No more complicated, clumsy remote controls. Enjoy more contrast and less jitter than a projector.

The Digital Eye Chart is ideal for offices using paperless patient records because it can use the same computer. Even with the purchase of a laptop or tablet, the Digital Eye Chart is less costly than purchasing a new projector.

**An affordable digital eye chart solution.**  
[www.DigitalEyeChart.com](http://www.DigitalEyeChart.com)  
 (208) 520 4292  
[sales@DigitalEyeChart.com](mailto:sales@DigitalEyeChart.com)

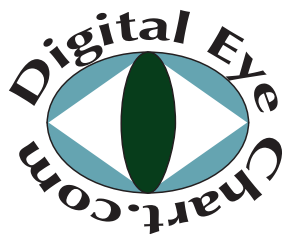

Many display options and chart types are included. Standard acuity and Red/Green charts are available with Letters, Numbers, Tumbling E's, Rolling C's, and Pictures, as well as, fixation, four dot, and astigmatism charts. Use the customizable slide show for promotions, education, and entertainment. If you need a different chart, let us know. Calibration Assistant (see bottom left) makes setup very easy.

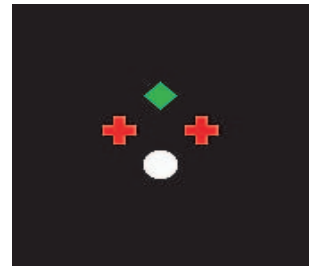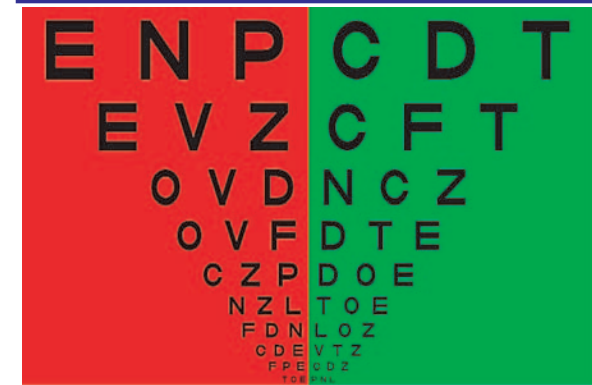


File Type Control

**E**

**Calibration Assistant**

Enter the distance from patient to wall chart. Click outside the entry box and the letter height will update. The 'E' used in calibration will be the size of the 20/100 acuity line on the wall chart. Adjust with the buttons above and measure the top 'E' on the wall chart.

Chart Distance	Calibration Letter Height			
	inches	in. decimal	in. fractions	mm
132		.96	0 - 31/32	24.4

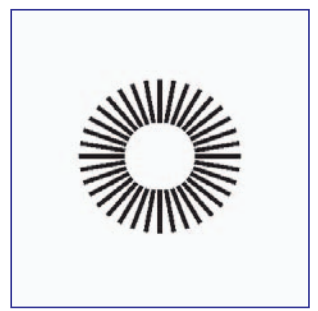

See [www.DigitalEyeChart.com](http://www.DigitalEyeChart.com) for details.